

Fortrex™
0.035" OTW PTA Balloon Catheter

Product code 135 cm catheter length	Product code 80 cm catheter length	Product code 40 cm catheter length	Balloon diameter (mm)	Balloon length (mm)	Nominal pressure (atm)	Rated burst pressure (atm)	Recommended introducer sheath (F)
A35HPV04020135	A35HPV04020080	A35HPV04020040	4 mm	20 mm	12	24	6
A35HPV04040135	A35HPV04040080	A35HPV04040040	4 mm	40 mm	12	24	6
A35HPV04080135	A35HPV04080080	A35HPV04080040	4 mm	80 mm	12	24	6
A35HPV04100135	A35HPV04100080	A35HPV04100040	4 mm	100 mm	12	24	6
A35HPV05020135	A35HPV05020080	A35HPV05020040	5 mm	20 mm	12	24	6
A35HPV05040135	A35HPV05040080	A35HPV05040040	5 mm	40 mm	12	24	6
A35HPV05080135	A35HPV05080080	A35HPV05080040	5 mm	80 mm	12	24	6
A35HPV05100135	A35HPV05100080	A35HPV05100040	5 mm	100 mm	12	24	6
A35HPV06020135	A35HPV06020080	A35HPV06020040	6 mm	20 mm	12	24	6
A35HPV06040135	A35HPV06040080	A35HPV06040040	6 mm	40 mm	12	24	6
A35HPV06080135	A35HPV06080080	A35HPV06080040	6 mm	80 mm	12	23	6
A35HPV06100135	A35HPV06100080	A35HPV06100040	6 mm	100 mm	12	23	6
A35HPV07020135	A35HPV07020080	A35HPV07020040	7 mm	20 mm	9	20	6
A35HPV07040135	A35HPV07040080	A35HPV07040040	7 mm	40 mm	9	20	6
A35HPV07080135	A35HPV07080080	A35HPV07080040	7 mm	80 mm	9	20	6
A35HPV07100135	A35HPV07100080	A35HPV07100040	7 mm	100 mm	9	20	6
A35HPV08040135	A35HPV08040080	A35HPV08040040	8 mm	40 mm	9	20	6
A35HPV08080135	A35HPV08080080	A35HPV08080040	8 mm	80 mm	9	19	6
A35HPV08100135	A35HPV08100080	A35HPV08100040	8 mm	100 mm	9	18	6
A35HPV09040135	A35HPV09040080	A35HPV09040040	9 mm	40 mm	9	18	7
A35HPV09080135	A35HPV09080080	A35HPV09080040	9 mm	80 mm	9	18	7
A35HPV10040135	A35HPV10040080	A35HPV10040040	10 mm	40 mm	8	16	7
A35HPV10080135	A35HPV10080080	A35HPV10080040	10 mm	80 mm	8	16	7
A35HPV12040135	A35HPV12040080	A35HPV12040040	12 mm	40 mm	8	14	7
A35HPV12080135	A35HPV12080080	A35HPV12080040	12 mm	80 mm	7	12	7



Indications, contraindications, warnings, and instructions for use can be found in the product labeling supplied with each device.

ENDURANCE UNDER PRESSURE

Fortrex™
0.035" OTW PTA Balloon Catheter



Medtronic

Medtronic Korea
17F, Glass Tower, #534
Teheran-ro, Gangnam-gu
Seoul, 06181, Korea
Tel: +82-2-3404-3150
Fax: +82-2-539-0961

UC201608108EE ©2016 Medtronic.
All rights reserved. Printed in Europe.

www.medtronic.com

Medtronic
Further Together

PREDICTABLE HIGH PRESSURE TREATMENT FOR AV ACCESS LESIONS

Providing you with an AV access solution at the value you need.

The high pressure Fortrex™ PTA balloon is the next-generation solution for deliverability, predictability, and procedural efficiency.

Deliverability Optimizing results.


- Low tip entry profile enables tight tracking on the wire
- Robust, flexible shaft design facilitates successful navigation in tortuous vessels

Predictability Even in challenging lesions.

- Balloon material and design allow for shape retention at higher pressure
- Focal pressure is exerted on the lesion for controlled, targeted treatment

Procedural efficiency Desired by physicians.

- Rapid deflation time contributes to reduced procedure length
- Balloon material and wall thickness permit reliable balloon rewrap into the sheath

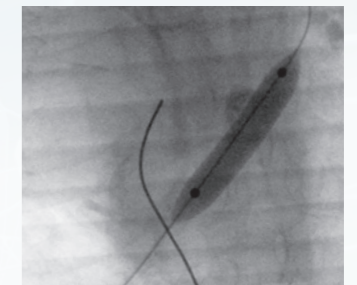

Fortrex balloon
Pressure is directed toward the lesion